



REHASENSE

## ICON Mobility System



## ICON 40

ALUMINUM LIGHTWEIGHT FOLDING  
WHEELCHAIR FOR ACTIVE USE.

**Icon 40** High performance, light weight, aluminium folding wheelchair with extensive set up options to meet the active Self Propelling user's needs.

# ICON 40

## Aluminum lightweight folding wheelchair for active use.

The **Icon 40** is part of the modular **Icon Mobility System** within the group of **Icon** wheelchairs with interchangeable parts. The **Icon 40** meets the wishes and requirements of a varied group of users. The many adjustment options offer the possibility to adapt the chair to the wishes of the therapist or the individual wishes of the user. The design has a sideward locking mechanism for the legrest to provide a clear front to avoid hooking and damaging of clothes, limbs or furniture.



REHASENSE



Legrests "swing-in swing-out, swing under", detachable. Flip up, angle adjustable footplates with heel band.



Armrests: with side guards, swing back for transfer, detachable and with height adjustable armpads.




Solid puncture proof rear wheel with grab ring. Easy to disconnect with the quick-release axle.



Angle adjustable front fork for perpendicular adjustment of the bearing house to the floor to optimize driving performance.



Slide up height adjustable anti-tip combined with step assist for easier tipping backwards.

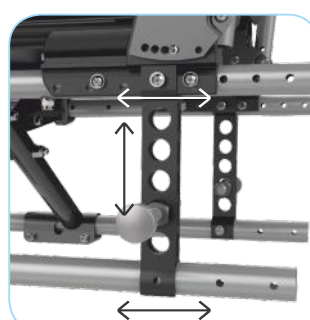

Multiple rear wheel positions for variation in seat depth and seat height.

### TECHNICAL DETAILS

Colour:  
Champagne-Black

Safe weight limit:  
150 kg / 190 kg (XL)

Weight:  
16,6 kg \*

Transportation weight:  
8,0 kg \*

Backrest height:  
42 cm

Seat height:  
48 cm

Total length:  
101 cm

Length without leg rests:  
76 cm


Height folded:  
92,5 cm

Width folded:  
34 cm

Seat widths:  
Standard: 40 - 50 cm  
XL: 53 - 61 cm

Seat depths:  
40 - 50 cm

\*) mid-size values

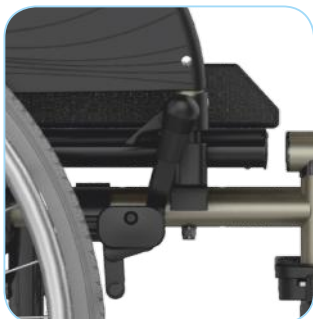

Standard hand brakes push to brake.



Velcro adjustable back for individual comfort.



Foldable cross brace for easier transportation.



Elevating legrests.  
Amputee legrests.  
Crutch holder.  
Height adjustable push handles.

### ACCESSORIES

The **Icon 40** is CE-marked and fulfils the requirements of EU Regulation 2017/745 (MDR) for medical devices and crash test of ISO 7176-19.

Copyright © Rehasense® 2021

Rehasense sp. z o.o.  
Sulejowska 45g  
97-300 Piotrków Trybunalski  
Poland